



PRODUCT INFORMATION

Product Name:

PROCEEDx Influenza A and Rotavirus Positive Sample

Catalogue Number:

VP-13-05

Product Description:

This product contains inactivated influenza A and rotavirus in storage buffer.

SPECIFICATIONS

Strains:

Influenza A/H1N1 Subtype

Simian rotavirus SA11, ICTV 00.060.0.03.001.00.201.001

Storage:

Store this preparation at 2 °C to 8 °C.

Shelf Life:

72 weeks from manufacturing date.

Volume:

1 mL

USE

Instructions For Use:

Use PROCEEDx product as per sample handling instructions specified by device or instrument manufacturer.

Intended Use:

The product is to be used for immunoassay instrument (device) verification/validation, manufacturing QMS, educational and R&D purposes.

NOT FOR USE IN DIAGNOSTIC PROCEDURES. FOR RESEARCH USE ONLY.

ASSISTANCE

If you have any questions regarding the production, testing or use of this antigen, please send them together with any relevant data to Microbix Technical Services, at customer.service@microbix.com or fax 905-361-8911. Your complete satisfaction with the performance of this product is important to us.

HAZARDS

We are aware of no specific biological hazards associated with this product. The preparation has been inactivated and should contain no infectious material however it should be handled as if potentially infectious. Generally accepted good laboratory practices appropriate to biological reagents should be employed. The buffer solution contains 0.02% sodium azide (CAS No. 26628-22-8) as a preservative.


SPECIAL HANDLING

No special handling required for this product.

QUALITY CONTROL INFORMATION



Test / Kit	Target	Result
TRU Influenza A & B Catalogue # 751230	Influenza A & B	Influenza A Positive Influenza B Negative
Binax Now Influenza A & B Catalogue # 416-000	Influenza A & B	Influenza A Positive Influenza B Negative
OSOM Influenza A & B Catalogue # 190	Influenza A & B	Influenza A Positive Influenza B Negative
BD Veritor Influenza A & B Catalogue # 256041	Influenza A & B	Influenza A Positive Influenza B Negative
SAS Rotavirus Catalogue # 69020	Rotavirus	Rotavirus Positive
Immuno Card STAT Rotavirus Catalogue # 750030	Rotavirus	Rotavirus Positive
Quickvue Chlamydia Catalogue # OB006	Chlamydia	Chlamydia Negative
Binax NOW RSV Catalogue # 430022	Respiratory Syncytial Virus	RSV Negative
SAS Adenovirus Catalogue # 067020	Adenovirus	Adenovirus Negative

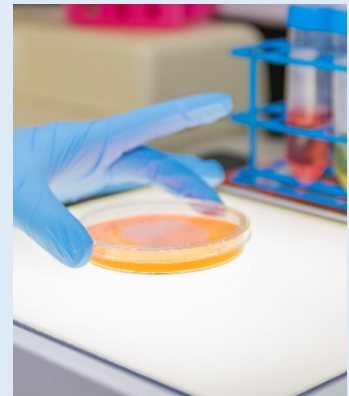

REF	Catalog Number
LOT	Lot Number
	Temperature Limitation
	Expiration Date
	Biological Risk